

## PET 18 PETROLEUM

## Analysis of Polywax 655 on Megabore BPX1

### KEYWORDS:


Polywax 655

Simulated distillation

Hydrocarbons

Part No.: 054800  
Phase: BPX1-0.1  
Column: 5m x 0.53mm ID  
Initial Temp: 40°C  
Rate: 15°C  
Final Temp: 420°C, 5 min  
Detector Temp: 440°C  
Carrier: He, 10mL/min  
Instrument: HP 6890  
Solvent: CS<sub>2</sub>

1. C20
2. C40
3. C60
4. C80
5. C100

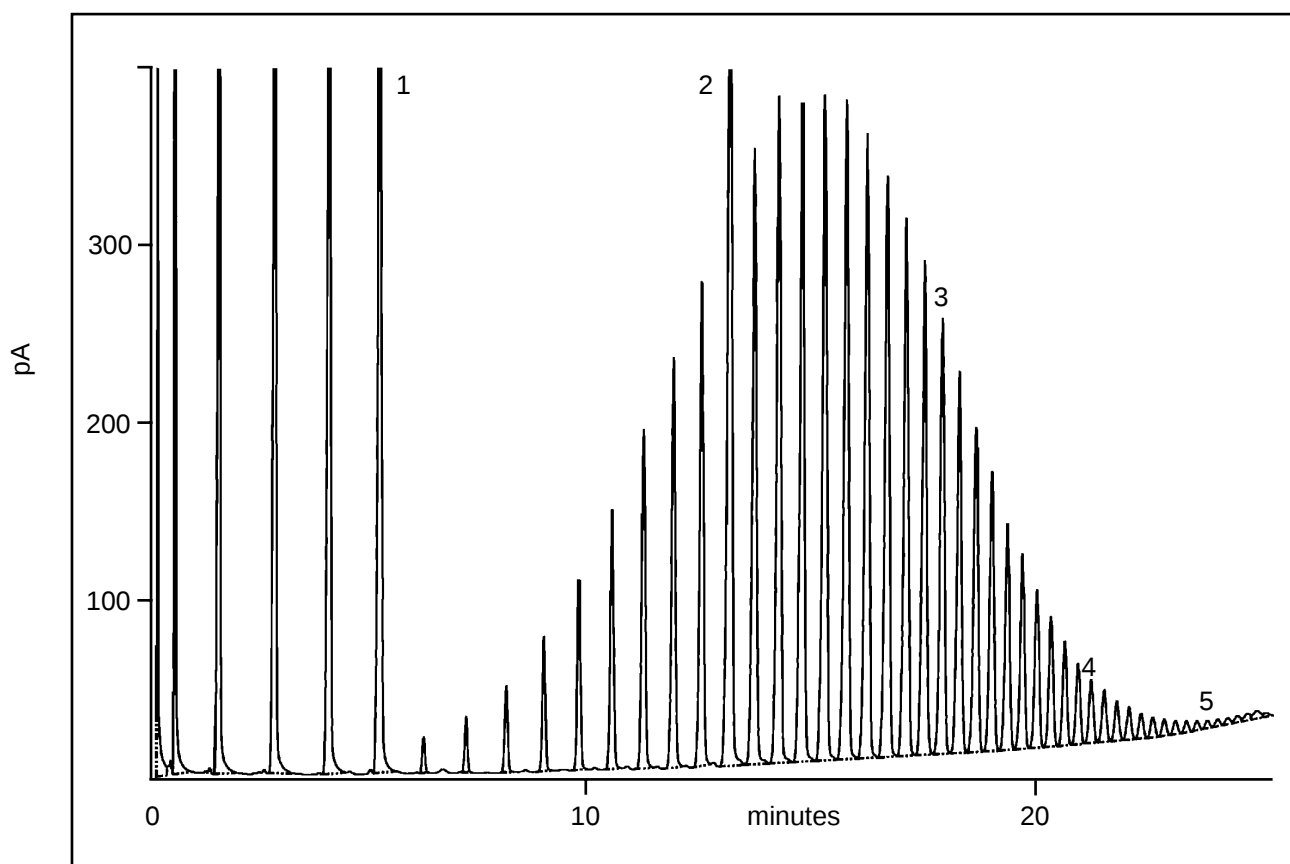

visit [www.sge.com](http://www.sge.com) or contact our technical Customer Support Team on: [techsupport@sge.com](mailto:techsupport@sge.com)

SGE International Pty. Ltd.  
email: [support@sge.com](mailto:support@sge.com)

SGE Incorporated (USA)  
email: [usa@sge.com](mailto:usa@sge.com)

SGE (France) Sarl  
email: [france@sge.com](mailto:france@sge.com)

SGE Deutschland GmbH  
email: [germany@sge.com](mailto:germany@sge.com)

SGE Japan Inc.  
email: [sgejpn@itjit.ne.jp](mailto:sgejpn@itjit.ne.jp)

SGE Australia Pty. Ltd.  
email: [sales@sge.com](mailto:sales@sge.com)

SGE Europe Ltd. (UK)  
email: [uk@sge.com](mailto:uk@sge.com)

SGE (Italia) Srl.  
email: [sge.italia@tin.it](mailto:sge.italia@tin.it)

SGE China Service Centre  
email: [sge@bj.col.com.cn](mailto:sge@bj.col.com.cn)

SGE India  
email: [sgeindia@bom4.vsnl.net.in](mailto:sgeindia@bom4.vsnl.net.in)

